


VENOSAN® 4000 RANGE

MEASURING GUIDE AND SIZE CHART

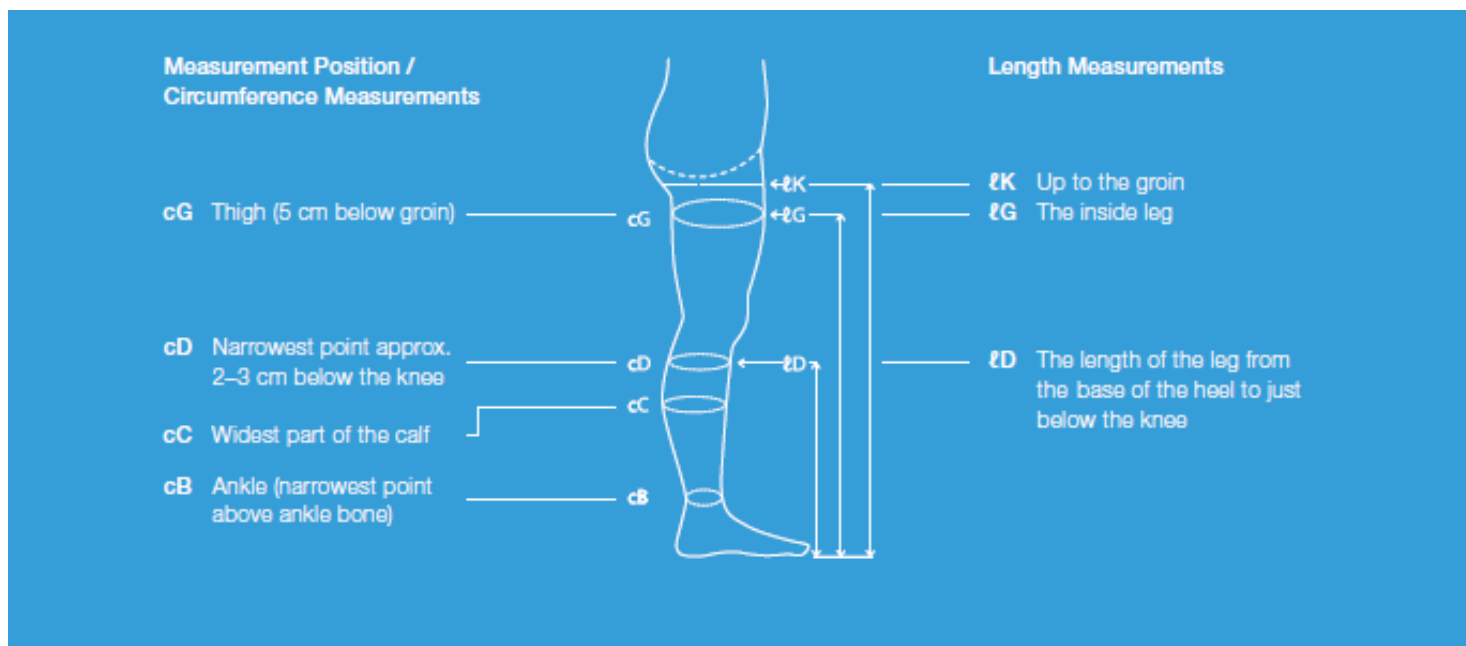

CIRCUMFERENCE MEASUREMENTS (CM)

	SMALL	MEDIUM	LARGE	X-LARGE
cB Ankle	20-23	23-26	26-29	29-31
cC Mid-calf	29-39	33-42	36-45	39-47
cD Below the knee	28-38	32-41	35-44	38-46
cG Thigh	48-63	54-67	58-72	63-74

LENGTH MEASUREMENTS (CM)

		short	long
AD/ADH	lD	35-40	37-42
AGH	lG	63-72	68-78
AGG/AT ATU Materna	lK	67-74	74-85

SIZE

It is important that we achieve the correct pressure gradient in accordance with the compression class prescribed by the doctor. The pressure exerted on the leg is determined by 2 factors:

- Physical law (LaPlace), which says that with constant tension the pressure is reduced with increased circumference and vice versa.
- With increased stretch the tension is increased and vice versa. One effect is therefore partially compensated through the other: Larger circumference = pressure loss due to circumference, but pressure gain due to increased tension of the material.

Thanks to more than 50 years experience we are familiar with these effects and can therefore offer a simple size range in accordance with our sizing chart, which guarantees medical efficacy.

For legs with a special shape or form we offer a made-to-measure service with extra fast delivery. Please ask for the special order form for made-to-measure stockings (REF 995220E).

LENGTH

VENOSAN® 4000 has a very high stretch in the length and adapts to almost every leg length without slipping. We have also proven through extensive testing that with a difference in length of +/- 3 cm, the pressure changes only by 3-5%.

Subsequently compare with the size chart for:

- Which size does the ankle circumference cB fall into (small, medium, large or x-large)?
- Which size do the other measurements fall into cC, cD and cG?

Custom made stockings may not be returned for credit. The manufacturer cannot be held responsible for incorrect measurements.